

instant accurate measurements

The new NSK Advanced Digital Apex Locator has been developed to accurately measure the length of canals in all conditions, including dry, wet and where blood is present.




- Advanced Digital Apex Locator
- User Friendly Interface
- Large Screen for Easy Monitoring
- Compact and Smart Design
- Audible Warning System



Large Display, Plus Audio Signals for Immediate Monitoring

The large iPex display allows immediate precise monitoring of the apical distance. User friendly graphics, plus audio signals, provides the accurate information required for the optimum file position as the instrument tip approaches the Apex. The audio volume is adjustable.

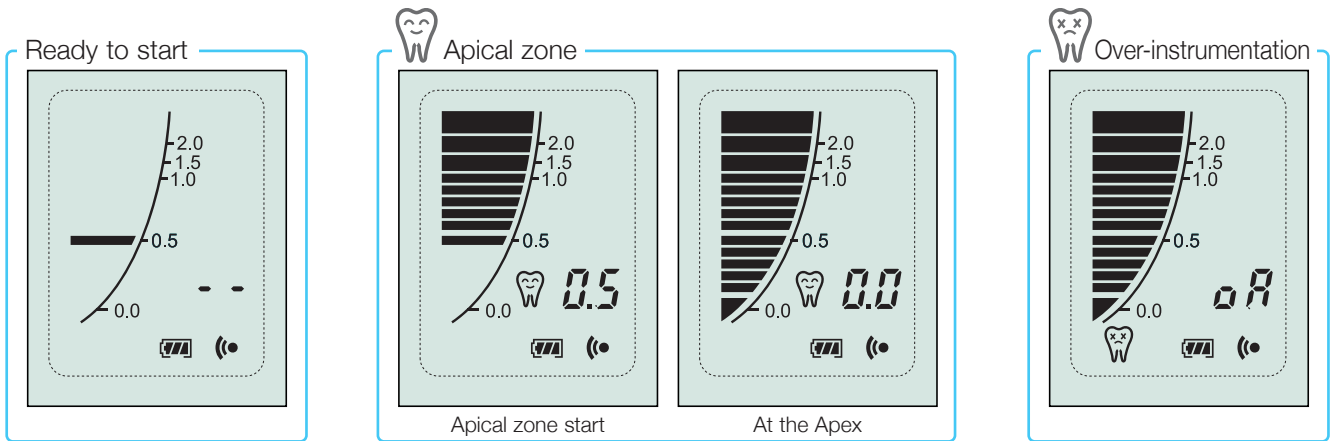

Accurate Digital Apex Locator

NSK's new digital apex locator insures great accuracy in all canal conditions and eliminates problems apparent in other Apex Locators which use analogue signals. The NSK iPex multi-frequency measuring software technology filters the irrelevant signals that can appear in curved canals or under other extraordinary circumstances. iPex also automatically select the best possible combinations of frequency to meet patient's canal condition.

- **Auto Power Shut-off**
iPex automatically shuts off after 10 minutes inactive.
- **Audible Warning System**
Audible signal informs clinicians about the progress of the file in the canal.
- **No Manual Calibration**
- **AAA Battery Operation**



iPex Complete Pack

| MODEL | ORDER CODE |
|-------|------------|
| iPex | Y141-251 |

Contents

- iPex unit
- Probe
- Extension cord
- Lip hook
- File clip
- AAA manganese dry batteries (4 pcs.)

Specifications

- Power Supply : for 1.5V AAA (LR03 or R03) batteries
- Continuous use : Approx. 50 hours
- Dimensions : W80 x D80 x H95 (mm)
- Weight : Unit 185g



File clip, lip hook and extension cord can be autoclaved at 135°C max. Two replacement file clips and lip hooks come as standard.



www.nsk-inc.com

NAKANISHI INC
700 Shimohinata
Kanuma, Tochigi
322-8666, Japan
TEL:+81 (0)289-64-3380
FAX:+81 (0)289-62-5636

NSK TOKYO OFFICE
6F MY Building
Ueno 1-13-3, Taito-ku
Tokyo 110-0005, Japan
TEL:+81 (0)3-3835-2892
FAX:+81 (0)3-3835-2856

NSK EUROPE GmbH EC REP
Westerbachstraße 58
D-60489, Frankfurt,
Germany
TEL:+49 (0)69 74 22 99 0
FAX:+49 (0)69 74 22 99 29

NSK FRANCE SAS
19 avenue de Villiers
75017 Paris,
France
TEL:+33 01 56 79 59 80
FAX:+33 01 56 79 59 81

NSK MIDDLE EAST
Dubai Airport Free Zone,
Room 2W-310,
3F West Wing, Dubai, UAE
TEL:+971-4-299-6644
FAX:+971-4-299-6664