

VARIOS 560 FROM NSK
MULTI-FUNCTION ULTRASONIC SYSTEM

NSK
Varios 560

compact and versatile

Now available with the Varios Optic handpiece, Varios 560 offers 3 power modes that can be selected at the touch of a button. Each mode corresponds to the various types of NSK tips. The choices are "G" (General), "E" (Endo) and "P" (Perio) mode. The selected mode is visually indicated via the LED display. The power button allows fine adjustment within each power mode setting.



Powerful Partners®

Lightweight and Compact

Weighs only 33g, Varios handpiece's slim and easy-to-handle body offers superior performance while reducing your fatigue even in long use. The vibration unit is built of solid titanium.

Circular Optic Illumination



The circular light emission of the VA-LUX-HP optic handpiece clearly illuminates the operational area

facilitating the possibility of faster treatment with reduced operator stress.

Automatic Tuning


When turned on, Varios handpieces automatically calibrate micromotors to the resistance of attached handpieces.

Wide Tip Range

Varios ultrasonic scaler tips are available for many clinical applications.



Both Optic and Non-optic are available for Varios 560

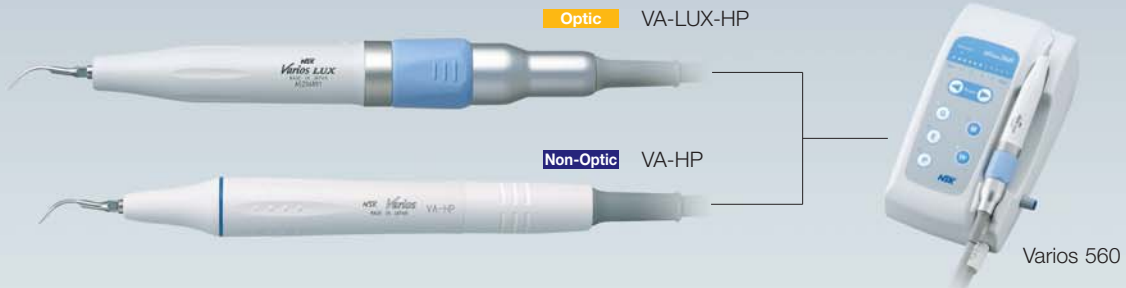

Varios 560 Complete Packs



Optic Varios 560 LUX Complete Pack

MODEL	ORDER CODE
VA560 LUX 120S	Y141-225
VA560 LUX 230S	Y141-228

Non-Optic Varios 560 Complete Pack

MODEL	ORDER CODE
VA560 120S	Y141-226
VA560 230S	Y141-229

Contents

- Control unit • Handpiece • Handpiece cord • 3 Scaling tips (G1, G4, G6)
- Water filter set • Water supply connector • Foot controller (FC-60)
- Tip replacement wrench with torque limiter • Tip holder • Power cord
- Tip cover S

Specifications

- Power Supply : AC 120V or 230V 50/60 Hz
- Frequency : 28~32 KHz
- Max. output : 8 W (G mode)
- Dimensions : W109 x D194 x H95 (mm)
- Unit weight : 1.2 kg

NSK Varios Ultrasonic Scaler Tip Book

NSK offers a wide range of ultrasonic tips for general scaling, endodontics, periodontics, restorative, crown loosening and plugging. The application for each tip is shown in the Varios TIP BOOK.



www.nsk-inc.com

NAKANISHI INC
700 Shimohinata
Kanuma, Tochigi
322-8666, Japan
TEL:+81 (0)289-64-3380
FAX:+81 (0)289-62-5636

NSK TOKYO OFFICE
6F MY Building
Ueno 1-13-3, Taito-ku
Tokyo 110-0005, Japan
TEL:+81 (0)3-3835-2892
FAX:+81 (0)3-3835-2856

NSK EUROPE GmbH EC REP
Westerbachstraße 58
D-60489, Frankfurt,
Germany
TEL:+49 (0)69 74 22 99 0
FAX:+49 (0)69 74 22 99 29

NSK FRANCE SAS
19 avenue de Villiers
75017 Paris,
France
TEL:+33 01 56 79 59 80
FAX:+33 01 56 79 59 81

NSK MIDDLE EAST
Dubai Airport Free Zone,
Room 2W-310,
3F West Wing, Dubai, UAE
TEL:+971-4-299-6644
FAX:+971-4-299-6664