

VARIOS 150 FROM NSK
BUILT-IN MODULE ULTRASONIC SCALERS

NSK
Varios 150

compact & reliable



Varios LUX Ultrasonic Scaler handpieces deliver optic illumination directly to the clinical treatment site.



Small Module for Integration Into Dental Equipment

Varios 150 is designed to be integrated into virtually any new or existing dental equipment. The scaler handpiece becomes included in the general handpiece console, and is controlled via the equipment foot control. Available in Optic and Non-Optic models.



Built-in Module



Optic

Non-Optic



Powerful Partners®

Varios 150 Standard Packs



Varios 150 LUX Complete Pack

Varios 150 LUX Complete Pack **Optic** Varios 150 Complete Pack **Non-Optic**

MODEL	ORDER CODE
VA150 LUX	Y141-010

MODEL	ORDER CODE
VA150	Y141-904

Contents

- Built in module • Handpiece • Handpiece cord • Handpiece holer • Power adjustment volume knob
- 3 Scaling tips (G1, G4, G6) • Tip replacement wrench with torque limiter • Tip holder


Specifications

- Power Supply : 24V AC • Frequency : 28-32 KHZ • Max. output : 8 W
- Dimensions : W57 x D73 x H33 (mm) • Unit weight : 114 g

Spare Parts

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
VA-LUX-HP Optic handpiece		E250-050	CR-1 Tip replacement wrench with torque limiter		Z221-070
VA-HP Non-optic handpiece		E177-050	E Tip Replacement Wrench		Z217-399
Varios 150 LUX Built-in Module		U387-012	Tip Holder		Z221-A080
VA 150 LUX-SC		E250-062	Tip Cover L		Z217-850
Varios 150 Built-in Module		U387-015	Tip Cover S		Z217-851
VA 150-SC		E177-062	VA Bulb (Pack of 3)		Y900-107

Varios 150 LUX PCB Connections

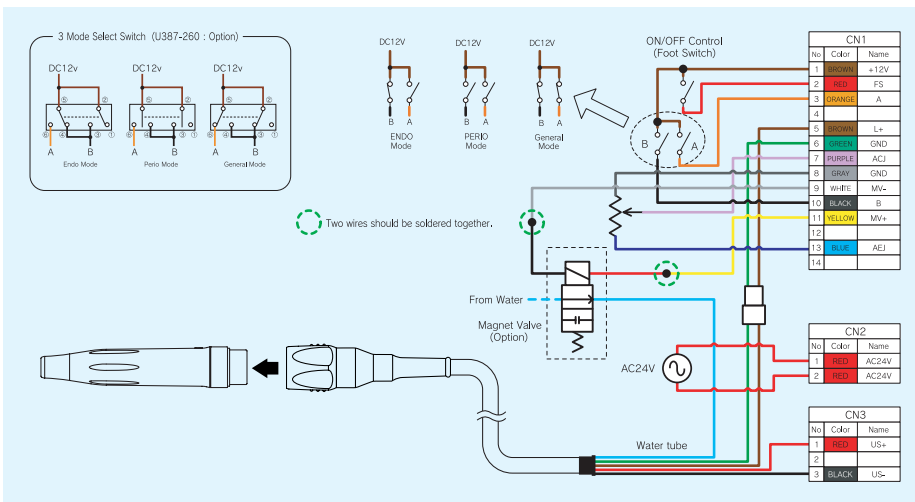

- The factory-set wiring allows only the Scaling (G-Mode).
- Connect the optional switch as shown above and choose the desired mode out of 3 (Perio, Endo, or Scaling) by this switch. (Optional switch : Part Number U387-260)

NSK Varios Ultrasonic Scaler Tip Book

NSK offers a wide range of ultrasonic tips for general scaling, endodontics, periodontics, restorative, crown loosening and plugging. The application for each tip is shown in the Varios TIP BOOK.



	www.nsk-inc.com					
	NAKANISHI INC 700 Shimohinata Kanuma, Tochigi 322-8666, Japan TEL:+81 (0)289-64-3380 FAX:+81 (0)289-62-5636	NSK TOKYO OFFICE 6F MY Building Ueno 1-13-3, Taito-ku Tokyo 110-0005, Japan TEL:+81 (0)3-3835-2892 FAX:+81 (0)3-3835-2856	NSK EUROPE GmbH <table border="1"> <tr> <td>EC</td> <td>REP</td> </tr> </table> Westerbachstraße 58 D-60489, Frankfurt, Germany TEL:+49 (0)69 74 22 99 0 FAX:+49 (0)69 74 22 99 29	EC	REP	NSK FRANCE SAS 19 avenue de Villiers 75017 Paris, France TEL:+33 01 56 79 59 80 FAX:+33 01 56 79 59 81
EC	REP					