

MACH-LITE XT FROM NSK
CELLULAR OPTIC HIGH SPEED HANDPIECES

NSK
MACH-LITE XT

stainless steel high performance handpieces



MACH-LITE XT Is New
Generation Cellular Optic
High-speed Handpiece
With Ceramic Ball Bearings.



Powerful Partners®

MACH-LITE XT

Optic

135°C
SSS


Mach-Lite XT is a new generation of cellular optic, triple water spray, stainless steel high speed handpiece. Ceramic ball bearings provide quieter performance at high speeds with less wear, resulting in greater reliability. They dramatically increased turbine life. The cellular glass optics provide 20% more illumination than fibre optics while illumination remains bright even after repeated autoclaving.

Clean Head System



Clean Head System extends turbine life, prevents debris from entering turbine to damage bearings. Cross-infection situations are also greatly reduced.

Ceramic Ball Bearings

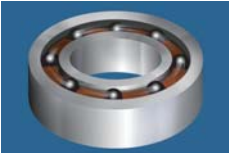

Ceramic Ball Bearings extends turbine life.

Cellular Optic Glass



Cellular Optic Glass is brighter and the most durable through repeated sterilizations.

Triple Water Spray



Triple Water Spray creates mist-like miniscule water droplets and sprays them evenly across the treatment field for maximum cooling effect.

Ultra Push Chuck



Ultra Push Chuck is NSK's original push button autochuck mechanism that firmly grips the bur during rotation.



MODEL ORDER CODE

MLXT-MU P621

- Miniature head

MLXT-SU P620

- Standard head

MLXT-TU P622

- Torque head

PTL-CL-4HV P401-051

- Light bulb included in coupling

PTL-CL-III P425-051

- Light bulb included in coupling
- With water volume adjuster

PTL-CL-FV P380-051

- Midwest 5 hole

Spare Parts

MODEL	PRODUCTS	ORDER CODE
NMC-MU03 Cartridge for Miniature head		P047
NMC-SU03 Cartridge for Standard head		P045
NMC-TU03 Cartridge for Torque head		P046
PTL O-ring Set (Pack of 5)		Y900-580

MODEL	PRODUCTS	ORDER CODE
MC-MU04 Head cap for Miniature head		P047-050
PTL II-MU04 Head cap for Standard head		P040-050
PTL II-TU04 Head cap for Torque head		P042-050
MC-SU05 for Miniature, Standard & Torque head		P040-765



www.nsk-inc.com

NAKANISHI INC
700 Shimohinata
Kanuma, Tochigi
322-8666, Japan
TEL:+81 (0)289-64-3380
FAX:+81 (0)289-62-5636

NSK TOKYO OFFICE
6F MY Building
Ueno 1-13-3, Taito-ku
Tokyo 110-0005, Japan
TEL:+81 (0)3-3835-2892
FAX:+81 (0)3-3835-2856

NSK EUROPE GmbH EC REP
Westerbachstraße 58
D-60489, Frankfurt,
Germany
TEL:+49 (0)69 74 22 99 0
FAX:+49 (0)69 74 22 99 29

NSK FRANCE SAS
19 avenue de Villiers
75017 Paris,
France
TEL:+33 01 56 79 59 80
FAX:+33 01 56 79 59 81

NSK MIDDLE EAST
Dubai Airport Free Zone,
Room 2W-310,
3F West Wing, Dubai, UAE
TEL:+971-4-299-6644
FAX:+971-4-299-6664