

200L / 300L / 400L FROM NSK

CELLULAR OPTIC HIGH SPEED HANDPIECES

NSK

200L/300L/400L

high precision - high performance
stainless steel handpieces

Connect Direct to
KaVo[®], Sirona[®], W&H[®], Bien-Air[®]
Coupling



Powerful Partners[®]



200L / 300L / 400L




NSK Patented Clean Head System



Anti-Contamination System virtually eliminates internal contamination of the handpiece by micro organisms. Clean Head system also extends turbine life by actively preventing debris from entering turbine during turbine rotation.

Ceramic Ball Bearings

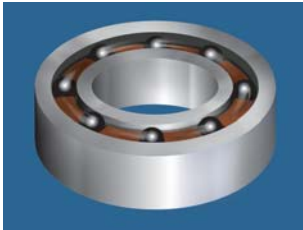

Quiet operation. Extended turbine life.

Cellular Glass Optics




20% brighter direct illumination of the work site. Will not deteriorate even after repeated sterilisation.

Triple Water Spray



Effectively delivers coolant spray direct to the work site.

NSK 200L / 300L / 400L is a series of superior high performance handpieces offering the long term durability of stainless steel bodies. Long term durability is further enhanced in the turbine with integrated Ceramic Bearings that deliver quiet concentric rotation with minimal wear, even after billions of revolutions. The overall result is increased reliability and intervals between scheduled technical service. Cellular glass optics provides optimum intraoral visibility through the transmission of 20% more illumination via the solid glass rod, as compared to conventional fibre optics which can deteriorate when individual fibres become redundant.

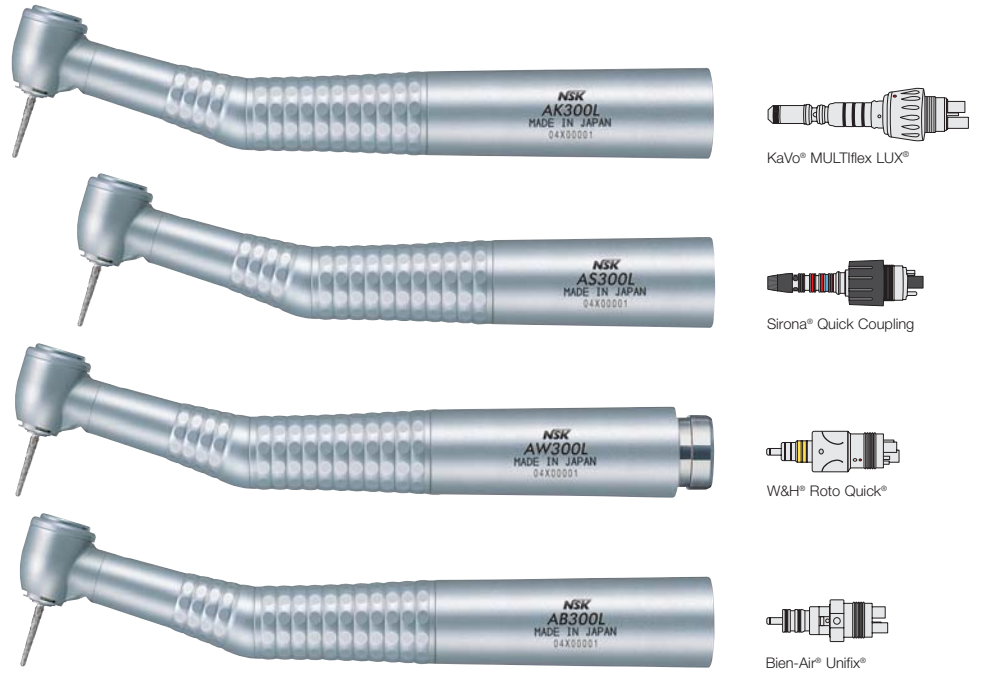

200L Miniature Head		300L Standard Head		400L Torque Head	
	13.1		14.65		14.65

	COUPLING TYPE	KaVo®	Sirona®	W&H®	Bien-Air®
200L Miniature Head	MODEL	AK200L	AS200L	AW200L	AB200L
	ORDER CODE	P671	P674	P688	P677
300L Standard Head	MODEL	AK300L	AS300L	AW300L	AB300L
	ORDER CODE	P672	P675	P689	P678
400L Torque Head	MODEL	AK400L	AS400L	AW400L	AB400L
	ORDER CODE	P673	P676	P690	P679

KaVo® and MULTiflex LUX® are registered trademarks of Kattenbach & Voigt GmbH & Co., Germany. Sirona® is a registered trademark of Sirona Dental Systems GmbH, Germany. W&H® and Roto Quick® are registered trademarks of W&H Dentalwerk Bürmoos GmbH, Austria. Bien-Air® and Unifix® are registered trademarks of Bien-Air Dental S.A., Switzerland.

Specifications are subject to change without notice.



All NSK handpieces are autoclavable up to 135°C.



www.nsk-inc.com

NAKANISHI INC.

700 Shimohinata
Kanuma, Tochigi
322-8666, Japan
TEL:+81 (0)289-64-3380
FAX:+81 (0)289-62-5636

NSK TOKYO Office

6F MY Building
Ueno 1-13-3, Taito-ku
Tokyo 110-0005, Japan
TEL:+81 (0)3-3335-2892
FAX:+81 (0)3-3335-2856

NSK EUROPE GmbH

Westerbachstraße 58
D-60489, Frankfurt
Germany
TEL:+49 (0)69 74 22 99 0
FAX:+49 (0)69 74 22 99 29

NSK FRANCE SAS

71 rue de Suresnes,
92380 Garches
France
TEL:+33 (0)1-4795-4517
FAX:+33 (0)1-4741-3516

NSK MIDDLE EAST

Dubai Airport Free Zone,
Room 2W-310,
3F West Wing, Dubai, UAE
TEL:+971-4-298-6644
FAX:+971-4-299-6664



Expect Perfection®