

Varios 350 Complete Packs



Specifications

- Power Supply : AC 120V or 230V 50/60 Hz
- Frequency : 28-32 KHz
- Max. output : 8 W (G mode)
- Dimensions : W80 x D140 x H32 (mm)
- Unit weight (Control unit, Handpiece) :
435g (350LUX), 393g (350)

Optic Varios 350 LUX Complete Pack

| MODEL | ORDER CODE |
|------------------|------------|
| VA350 LUX (120V) | Y140-960 |
| VA350 LUX (230V) | Y140-962 |

Non-Optic Varios 350 Complete Pack

| MODEL | ORDER CODE |
|--------------|------------|
| VA350 (120V) | Y140-952 |
| VA350 (230V) | Y140-954 |

Contents

- Control unit (with handpiece cord) • AC adaptor • Handpiece • Handpiece holder • Water supply tube with water filter • Water supply connector
- Foot pedal (FC-61) • 3 Scaling tips (G1, G4, G6) • Tip replacement wrench with torque limiter • Tip holder

Spare Parts

| MODEL | PRODUCTS | ORDER CODE | MODEL | PRODUCTS | ORDER CODE |
|---|----------|--------------|---|----------|------------|
| VA-LUX-HP Optic handpiece | | E250-050 | Tube Ring (Pack of 2) | | Y900-077 |
| VA-HP Non-optic handpiece | | E177-050 | CR-1 Tip replacement wrench with torque limiter | | Z221-070 |
| Varios 350 LUX Unit | | U449 | E Tip Replacement Wrench | | Z217-399 |
| Varios 350 Unit | | U423 | Tip Holder | | Z221-A080 |
| VA350 AC Adaptor (120V) VA350 AC Adaptor (230V) For Varios350 LUX & 350 | | U429 U430 | Tip Cover L | | Z217-850 |
| FC-61 Foot Control for Varios 350 LUX & 350 (4m) | | Z300-001 | Tip Cover S | | Z217-851 |
| Water Supply Tube | | U387-580 | Varios 350 Stand | | Z254 |
| Water Supply Connector | | U387-030 | Varios 350 Mounting Bracket | | Z259 |
| Water Filter Set (Irrigation tube, water filter, water filter case) | | U387-040 | Unit Holder For Varios 350 | | Z253 |
| Water Filter Case | | Y900-069 | VA Bulb (Pack of 3) | | Y900-107 |
| Water Filter | | U387-042 | | | |

VARIOS 350 LUX & VARIOS 350 FROM NSK MULTI-FUNCTION ULTRA SONIC SCALERS

NSK
Varios 350

compact & easy



NSK www.nsk-inc.com

NAKANISHI INC
700 Shimohinata
Kanuma, Tochigi
322-8666, Japan
TEL:+81 (0)289-64-3380
FAX:+81 (0)289-62-5636

NSK TOKYO OFFICE
6F MY Building
Ueno 1-13-3, Taito-ku
Tokyo 110-0005, Japan
TEL:+81 (0)3-3835-2892
FAX:+81 (0)3-3835-2856

NSK EUROPE GmbH
Westerbachstraße 58
D-60489, Frankfurt,
Germany
TEL:+49 (0)69 74 22 99 0
FAX:+49 (0)69 74 22 99 29

NSK FRANCE SAS
19 avenue de Villiers
75017 Paris,
France
TEL:+33 01 56 79 59 80
FAX:+33 01 56 79 59 81

NSK MIDDLE EAST
Dubai Airport Free Zone,
Room 2W-310,
3F West Wing, Dubai, UAE
TEL:+971-4-299-6644
FAX:+971-4-299-6664

REF: PR-D207E Ver.5 06.05.05.Y

Specifications are subject to change without notice.

Powerful Partners®

Compact Size Control Unit & Lightest Handpiece in the World.

Varios 350 has been developed especially to offer a compact size portable control unit, allowing easy installation into any dental unit. It is simple to operate and suitable for a variety of clinical procedures - scaling, perio and endo. Available in optic and non-optic models.



Simple design and easy procedure selection control.

One dial covers full procedure range (Scaling, Perio and Endo).



At a turn of the dial, any one of 3 power modes can be selected, corresponding to the types of tips. The choices are "G" (General), "E" (Endo) and "P" (Perio) mode.

The chosen mode can be visually checked via the colour coded indicator on the module fascia. The power dial allows easy power adjustment. Additionally, by selecting the "P" mode, the absolute minimum stable power is activated, enabling the most efficient perio treatment.

Allows Simple Installation Into All Dental Units.

Varios 350 Units Installation Examples



Right side of the dental mobile cart.

Behind the rear of the dental mobile cart.

Under the dental mobile cart.

On the dental console arm.

Ultrasonic Scalers Offering "One Does All" Versatility

In combination with an extensive range of tips, Varios Ultrasonic Scaler series introduce unprecedented "One handpiece does all" versatility for your every clinical application.



Optic VA-LUX-HP Optic Handpiece



Non-Optic VA-HP Non-Optic Handpiece


Lightweight and Compact

The small, easy-to-hold handpiece weighs just 33g (VA-HP), the lightest ultrasonic scaler handpiece in the world. Treatment is easy, and hand fatigue after a prolonged use is minimized.

Circular Optic Illumination



The circular light emission of the VA-LUX-HP optic handpiece clearly illuminates the operational area facilitating the possibility of faster treatment with reduced operator stress.

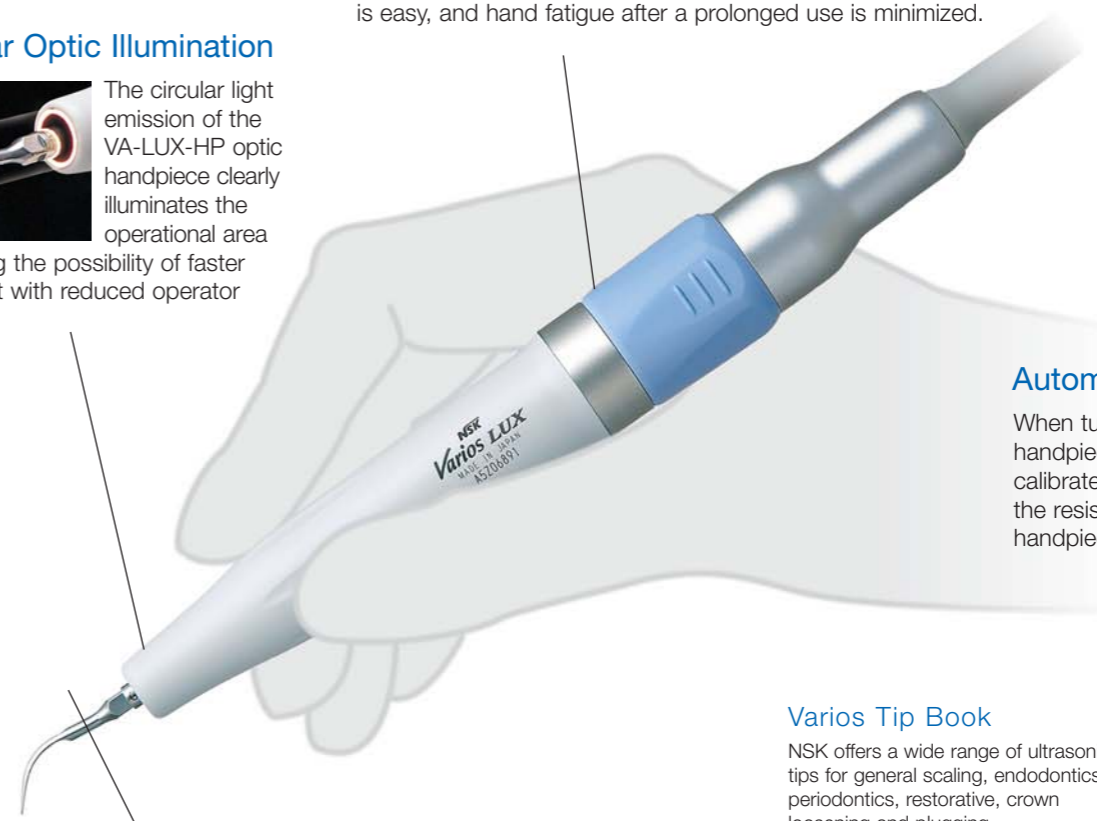

Automatic Tuning

When turned on, Varios handpieces automatically calibrate micromotors to the resistance of attached handpieces.

Wide Tip Range

Varios ultrasonic scaler tips are available for many clinical applications.



Varios Tip Book

NSK offers a wide range of ultrasonic tips for general scaling, endodontics, periodontics, restorative, crown loosening and plugging. The application for each tip is shown in the Varios TIP BOOK.

