






## PRESTO AQUA II Spare Parts

MODEL	PRODUCTS	ORDER CODE
PRESTO AQUA II Unit		Z308-015
PRA-HP PRESTO AQUA Handpiece		T816
Hose for PRESTO AQUA II		Z308-050
QD-J B3 Coupling		T354-051
AFC-40 Foot Control		Z127
PR-AQ03 PRESTO AQUA Cartridge		T816-016
Chuck		T815-019
Secondary Hose for PRESTO (2m)		Z101-240
Water Filter Set (water tube x 2, water filter, water filter case)		U387-040

MODEL	PRODUCTS	ORDER CODE
Water Filter with Case		Y900-069
Water Filter		U387-042
Tube Ring (includes 2)		Y900-077
PRA II Bottle Set		Z308-040
PRA II Bottle Lid		Z308-041
PRA II Bottle		Z308-601
PRA II Bottle O-ring Set (pack of 2)		Y900-327
Bracket Set		Z308-060
Chuck Wrench		C034-751

## PRESTO II Spare Parts

MODEL	PRODUCTS	ORDER CODE
PRESTO II Unit		Z307-015
PR-HP PRESTO Handpiece (includes 2m hose)		T815
AFC-40 Foot Control		Z127
Handpiece Stand		Z010-134
PR-03 PRESTO Cartridge		T815-016

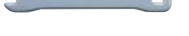MODEL	PRODUCTS	ORDER CODE
Chuck		T815-019
Hose for PRESTO (2m)		T815-220
Secondary Hose for PRESTO (2m)		Z101-240
LAB Spanner Wrench		H005-751
Chuck Wrench		C034-751

# NSK

## LUBRICATION FREE AIR TURBINE

## PRESTO AQUA II PRESTO II

# precision for professionals



Ideal for preventing  
ceramic fractures

PRESTO AQUA II



PRESTO II



Powerful Partners®

**NAKANISHI INC.**  
www.nsk-inc.com  
700 Shimohinata Kanuma-shi  
Tochigi 322-8666,  
Japan

TEL : +81 (0)289-64-3380  
FAX : +81 (0)289-62-5636

**NSK Europe GmbH**  
www.nsk-europe.de  
Westerbachstraße 58  
D-60489 Frankfurt,  
Germany

TEL : +49 (0)69 74 22 99 0  
FAX : +49 (0)69 74 22 99 29

**NSK United Kingdom Ltd.**  
www.nsk-inc.com  
Gateway1000, Arlington Business Park,  
Whittle Way, Stevenage,  
Hertfordshire, SG1 2FP, UK

TEL : +44 800 634 1909  
FAX : +44 800 634 1910

**NSK America - Canadian Branch**  
www.nsk-inc.com  
1054 Centre St. - Suite 283  
Thornhill, Ontario, Canada  
L4J 8E5

TEL/FAX : +1 905 780 6757  
Toll Free : +1 866 780 6757

**NSK Oceania Pty Ltd.**  
www.nsk-inc.com  
Unit 22, 198-222 Young St.  
Waterloo, Sydney,  
NSW 2017 Australia

TEL : +61 (0)2 8306 3000  
FAX : +61 (0)2 8306 3030

**NSK Middle East**  
www.nsk-inc.com  
Dubai Airport Free Zone,  
Room 2W-310,  
3F West Wing, Dubai, UAE

TEL : +971 (0)4-299-6644  
FAX : +971 (0)4-299-6664

# PRESTO AQUA II

Lubrication free laboratory air turbine handpiece with water spray

The PRESTO AQUA II comes with a coolant water spray that eliminates heat buildup on target substrate to prevent micro-cracks and decrease dispersal of debris and chips. An efficient tool for all ceramic work, it is particularly well suited to zirconium-based ceramics. The PRESTO AQUA II never requires lubrication, a feature that not only saves materials from contamination by oil but also contributes to a tidy work environment.



The PRESTO AQUA II delivers precision and ease to ceramic cutting by eliminating heat with a coolant water spray.



### Two Water Spray Options

The water spray can be switched between droplets or fine mist spray.



**Special Dust-proof Mechanism**  
Our original dust-proof mechanism prevents the infiltration of chips into the handpiece bearings to extend cartridge life.



**Easy-release Bur Chuck**  
The bur can be removed with ease with the twist of a ring.

## PRESTO AQUA II Complete Pack

MODEL	ORDER CODE
PRESTO AQUA II Complete pack	Y150-023

- Contents
- PRESTO AQUA II unit
  - PRESTO AQUA handpiece
  - QD-J B3 coupling
  - Hose for PRESTO AQUA II
  - Foot control (AFC-40)
  - Secondary hose for PRESTO


- Specifications
- 320,000min<sup>-1</sup>
  - Coolant water spray
  - Bottle/tap water switching function
  - Special dust-proof mechanism
  - No lubrication required

## Unique, easy bottle removal function



The water bottle (250mL) is easily removed with a release button. Operators can use the original water bottle or press the button to switch to tap water (direct water system).

## PRESTO AQUA handpiece

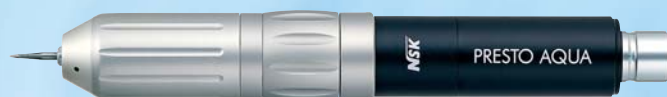

MODEL	ORDER CODE
PRA-HP	T816



MODEL	ORDER CODE
QD-J B3	T354-051

# PRESTO II

Lubrication free laboratory air turbine handpiece

This user-friendly air turbine handpiece is lightweight and compact, measuring just 130mm in length. Its ergonomically designed grip assures user comfort and easy operation even during demanding procedures.

## PRESTO II Complete Pack

MODEL	ORDER CODE
PRESTO II Complete pack	Y150-022

- Contents
- PRESTO II unit
  - PRESTO handpiece
  - Handpiece stand
  - Foot control (AFC-40)
  - Hose for PRESTO
  - Secondary hose for PRESTO

- Specifications
- 320,000min<sup>-1</sup>
  - Special dust-proof mechanism
  - No lubrication required
  - Water spray not employed

## PRESTO handpiece



MODEL	ORDER CODE
PR-HP	T815

