

## Varios 350 Complete Packs



### Specifications

- Power Supply : AC 120V or 230V 50/60 Hz
- Frequency : 28-32 KHz
- Max. output : 8 W (G mode)
- Dimensions : W80 x D125 x H32 (mm)
- Unit weight (Control unit, Handpiece) : 440g (350LUX), 393g (350)

### Optic Varios 350 LUX Complete Pack

MODEL	ORDER CODE
VA350 LUX (120V)	Y140-960
VA350 LUX (230V)	Y140-962

### Non-Optic Varios 350 Complete Pack

MODEL	ORDER CODE
VA350 (120V)	Y140-952
VA350 (230V)	Y140-954

### Contents

- Control Unit with Handpiece Cord • AC Adaptor • Handpiece • Handpiece Holder • Water Filter Set • Water Supply Connector
- Foot Control (FC-61) • 3 Scaling Tips (G1, G4, G6) • Tip Replacement Wrench with Torque Limiter • Tip Holder • Tip Cover S

## Spare Parts

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
VA-LUX-HP Optic Handpiece		E250-050	Tube Ring (Pack of 2)		Y900-077
VA-HP Non-Optic Handpiece		E177-050	CR-1 Tip Replacement Wrench with Torque Limiter		Z221-070
Varios 350 LUX Unit		U449	E Tip Replacement Wrench		Z217-399
Varios 350 Unit		U423	Tip Holder		Z221-080
VA350 AC Adaptor for Varios350 LUX & 350		U429-001(120V) U430-001(230V)	Tip Cover L		Z217-850
FC-61 Foot Control (4m) for Varios 350 LUX & 350		Z300-001	Tip Cover S		Z217-851
Water Supply Tube (2m)		U387-580	Varios 350 Stand		Z254
Water Supply Connector		U387-030	Varios 350 Mounting Bracket		Z259
Water Filter Set (2 water supply tube, water filter, water filter case)		U387-040	Unit Holder for Varios 350		Z253
Water Filter Case		Y900-069	VA Bulb (Pack of 3)		Y900-107
Water Filter		U387-042	Handpiece Holder for Varios 350		Z250-001
Middle Sheath		E250-075			

# NSK

## MULTI-FUNCTION ULTRA SONIC SCALERS

# Varios 350

# compact & easy installation



**NAKANISHI INC.**  
www.nsk-inc.com  
700 Shimohinata Kanuma-shi  
Tochigi 322-8666,  
Japan

TEL : +81 (0)289-64-3380  
FAX : +81 (0)289-62-5636

**NSK Europe GmbH**   
www.nsk-europe.de  
Westerbachstraße 58  
D-60489 Frankfurt,  
Germany

TEL : +49 (0)69 74 22 99 0  
FAX : +49 (0)69 74 22 99 29

**NSK United Kingdom Ltd.**  
www.nsk-uk.com  
Office 5, Gateway1000,  
Arlington Business Park,  
Whittle Way, Stevenage, SG1 2FP, UK

TEL : +44 (0)1438 310670  
FAX : +44 (0)1438 310671

**NSK Oceania Pty Ltd.**  
www.nsk-inc.com  
Unit 22, 198-222 Young St.,  
Waterloo, Sydney,  
NSW 2017 Australia

TEL : +61 (0)2 8306 3000  
FAX : +61 (0)2 8306 3030

**NSK America - Canadian Branch**  
www.nsk-inc.com  
1054 Centre St. - Suite 283  
Thornhill, Ontario, Canada  
L4J 8E5

TEL/FAX : +1 905 780 6757  
Toll Free : +1 866 780 6757

**NSK Middle East**  
www.nsk-inc.com  
Dubai Airport Free Zone,  
Room 2W-310,  
3F West Wing, Dubai, UAE

TEL : +971 (0)4-299-6644  
FAX : +971 (0)4-299-6664



## Compact Size Control Unit & Lightest Handpiece

Varios 350 has been developed especially to offer a compact size portable control unit, allowing easy installation into any dental unit. It is simple to operate and suitable for a variety of clinical procedures - Scaling, Perio and Endo. Available in Optic and Non-Optic models.

### One dial covers full procedure range (Scaling, Perio, and Endo).

At a turn of the dial, any one of 3 power modes can be selected, corresponding to the types of tips. The choices are "G" (General), "E" (Endo) and "P" (Perio) mode. The chosen mode can be visually checked via the colour coded indicator on the module facia. The power dial allows easy power adjustment. Additionally, by selecting the "P" mode, the absolute minimum stable power is activated, enabling the most efficient Perio treatment.



## Ultrasonic Scalers Offering "One Does All" Versatility

In combination with an extensive range of tips, Varios Ultrasonic Scaler series introduce unprecedented "One handpiece does all" versatility for your every clinical application.



**Optic** VA-LUX-HP Optic Handpiece



**Non-Optic** VA-HP Non-Optic Handpiece


### Lightweight and Compact

The small, easy-to-hold handpiece weighs just 33g (VA-HP), the lightest ultrasonic scaler handpiece in the world. Treatment is easy, and hand fatigue after a prolonged use is minimized.

### Circular Optic Illumination



The circular light emission of the VA-LUX-HP optic handpiece clearly illuminates the operational area facilitating the possibility of faster treatment with reduced operator stress.

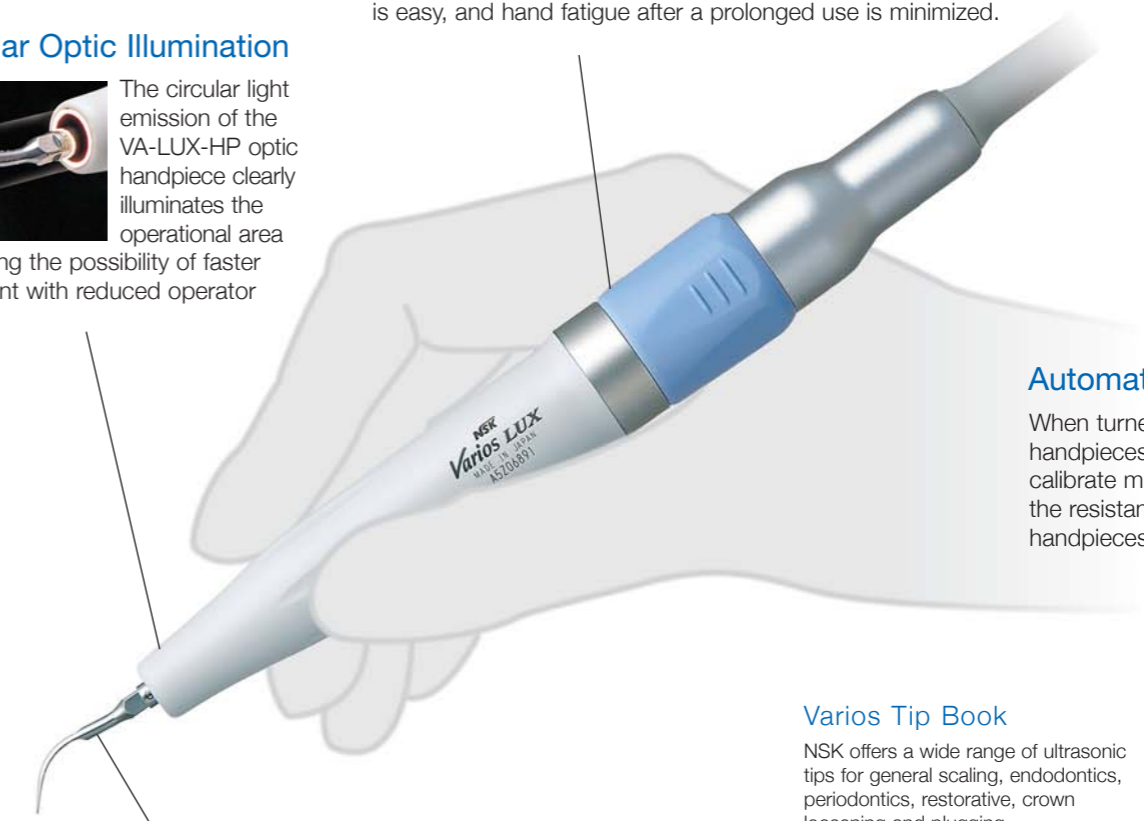

### Automatic Tuning

When turned on, Varios handpieces automatically calibrate micromotors to the resistance of attached handpieces.

### Varios Tip Book

NSK offers a wide range of ultrasonic tips for general scaling, endodontics, periodontics, restorative, crown loosening and plugging. The application for each tip is shown in the Varios TIP BOOK.



### Wide Tip Range

Varios ultrasonic scaler tips are available for many clinical applications.



### Allows Simple Installation Into All Dental Units. Varios 350 Units Installation Examples



Right side of the dental mobile cart.



Behind the rear of the dental mobile cart.



Under the dental mobile cart.



On the dental console arm.