

compact & reliable

Small Module for Integration Into Dental Equipment

Varios 150 is designed to be integrated into virtually any new or existing dental equipment. The scaler handpiece becomes included in the general handpiece console, and is controlled via the equipment foot control. Available in Optic and Non-Optic models.

Varios LUX Ultrasonic Scaler handpieces deliver optic illumination directly to the clinical treatment site.



Built-in Module



VA-LUX-HP

Optic




VA-HP

Non-Optic



Varios 150 Standard Packs

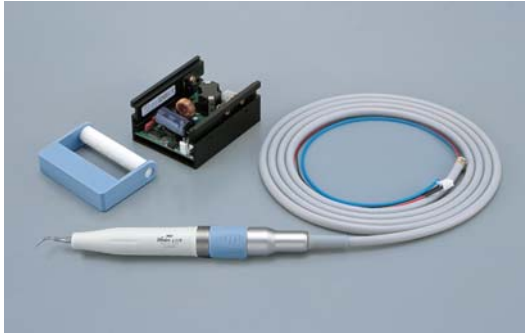

Varios 150 LUX Complete Pack **Optic** Varios 150 Complete Pack **Non-Optic**

MODEL	ORDER CODE
VA150 LUX	Y141-010

MODEL	ORDER CODE
VA150	Y140-904

Contents

- Built in module • Handpiece • Handpiece cord • Handpiece holder • Power adjustment volume knob
- 3 Scaling tips (G1, G4, G6) • Tip replacement wrench with torque limiter • Tip Holder


Specifications

- Power Supply : 24V AC • Frequency : 28-32 KHz • Max. output : 8 W
- Dimensions : W57 x D73 x H33 (mm) • Unit weight : 114 g

Spare Parts

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
VA-LUX-HP Optic handpiece		E250-050	CR-1 Tip replacement wrench with torque limiter		Z221-070
VA-HP Non-optic handpiece		E177-050	E Tip Replacement Wrench		Z217-399
Varios 150 LUX Built-in Module		U387-012	Tip Holder		Z221-080
VA 150 LUX-SC		E250-062	Tip Cover L		Z217-850
Varios 150 Built-in Module		U387-015	Tip Cover S		Z217-851
VA 150-SC		E177-062	VA Bulb (Pack of 3)		Y900-107

Varios 150 LUX PCB Connections

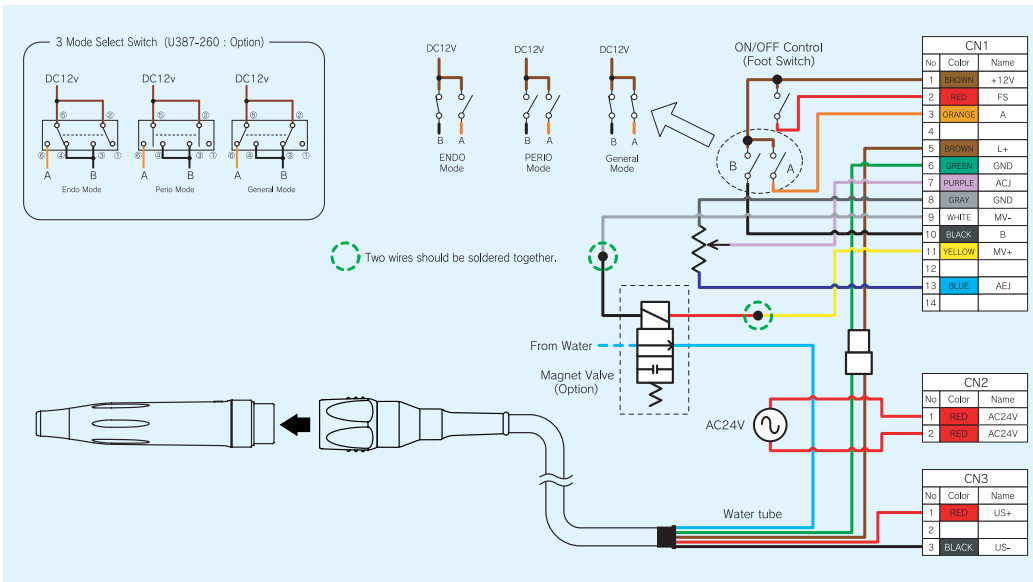

- The factory-set wiring allows only the Scaling (G-Mode).
- Connect the optional switch as shown above and choose the desired mode out of 3 (Perio, Endo, or Scaling) by this switch. (Optional switch : Part Number U387-260)

NSK Varios Ultrasonic Scaler Tip Book

NSK offers a wide range of ultrasonic tips for general scaling, endodontics, periodontics, restorative, crown loosening and plugging. The application for each tip is shown in the Varios TIP BOOK.



NAKANISHI INC.
www.nsk-inc.com
700 Shimohinata Kanuma-shi
Tochigi 322-8666,
Japan
TEL : +81 (0)289-64-3380
FAX : +81 (0)289-62-5636

NSK Europe GmbH EC REP
www.nsk-europe.de
Westerbachstraße 58
D-60489 Frankfurt,
Germany
TEL : +49 (0)69 74 22 99 0
FAX : +49 (0)69 74 22 99 29

NSK United Kingdom Ltd
www.nsk-uk.com
Office 5, Gateway1000,
Arlington Business Park, Whittle Way,
Stevenage, SG1 2FP, UK
TEL : +44 (0)1438 310670
FAX : +44 (0)1438 310671

NSK Oceania Pty Ltd
www.nsk-inc.com
Unit 22, 198-222 Young St.
Waterloo, Sydney,
NSW 2017, Australia
TEL : +61 (0)2 8306 3000
FAX : +61 (0)2 8306 3030

NSK America - Canadian Branch
www.nsk-inc.com
1054 Centre St, Suite 283
Thornhill, Ontario,
Canada L4J 8E5
TEL/FAX : +1 905 780 6757
Toll Free : +1 866 780 6757

NSK Middle East
www.nsk-inc.com
Dubai Airport Free Zone,
East Wing, Room 6EA-701,
PO Box 54316 Dubai, UAE
TEL : +971 (0)4-299-6644
FAX : +971 (0)4-299-6664