

# PANA-MAX **Non-Optic**

- Clean Head System
- Ultra push chuck
- Single Spray



MODEL **PAX-SU M4** ORDER CODE **P560-001**

- Standard Head
- Midwest 4 hole



MODEL **PAX-TU M4** ORDER CODE **P561-001**

- Torque Head
- Midwest 4 hole



MODEL **PAX-SU B2** ORDER CODE **P562-001**

- Standard Head
- Borden 2 hole



MODEL **PAX-TU B2** ORDER CODE **P563-001**

- Torque Head
- Borden 2 hole

## Spare Parts

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
<b>PAX-SU03</b>	Cartridge for standard head, includes O-ring & wave washer	Y900-212	<b>PAX-SU04</b>	Head cap for standard head	P074-050	<b>PAX-SU05</b>	Head cap wrench for standard head & torque head	P074-765
<b>PAX-TU03</b>	Cartridge for torque head, includes O-ring & wave washer	Y900-213	<b>PAX-TU04</b>	Head cap for torque head	P075-050			



Washable in the thermoisinfector.



Autoclavable at 135°C max.

**NAKANISHI INC.**  
www.nsk-inc.com  
700 Shimohinata Kanuma-shi  
Tochigi 322-8666,  
Japan

TEL : +81 (0)289-64-3380  
FAX : +81 (0)289-62-5636

**NSK Europe GmbH**  
www.nsk-europe.de  
Westerbachstraße 58  
D-60489 Frankfurt,  
Germany

TEL : +49 (0)69 74 22 99 0  
FAX : +49 (0)69 74 22 99 29

**NSK United Kingdom Ltd.**  
www.nsk-inc.com  
Gateway1000, Arlington Business Park,  
Whittle Way, Stevenage,  
Hertfordshire, SG1 2FP, UK

TEL : +44 800 634 1909  
FAX : +44 800 634 1910

**NSK Oceania Pty Ltd.**  
www.nsk-inc.com  
Unit 22, 198-222 Young St.  
Waterloo, Sydney,  
NSW 2017 Australia

TEL : +61 (0)2 8306 3000  
FAX : +61 (0)2 8306 3030

**NSK America - Canadian Branch**  
www.nsk-inc.com  
1054 Centre St. - Suite 283  
Thornhill, Ontario, Canada  
L4J 8E5

TEL/FAX : +1 905 780 6757  
Toll Free : +1 866 780 6757

**NSK Middle East**  
www.nsk-inc.com  
Dubai Airport Free Zone,  
Room 2W-310,  
3F West Wing, Dubai, UAE

TEL : +971 (0)4-299-6644  
FAX : +971 (0)4-299-6664

# NSK

## HIGH SPEED HANDPIECES

# PANA-MAX

# stunning power and durability



Powerful Partners®

# Brand-New PANA-MAX

NSK's Best-selling High Speed Handpiece Gets Even Better



## 1 Highest Torque in its Class

With robust 16.5-watt torque, the new PANA-MAX high speed handpiece delivers the best drilling performance within this category. Substantially improved low noise level, together with newly designed Arc-Line Body, promises ease of use and comfort for both patients and clinicians.

## 2 Clean Head System

The patented NSK Clean Head System is a special mechanism designed to prevent air retraction, thereby preventing the entry of oral fluids, viruses, and other foreign material into the handpiece head. An additional advantage of the Clean Head System is to prolong the life of the turbine bearings by blocking the entry of micro-contaminants into the turbine head.

## 3 NSK Ultra Push Chuck

NSK's patented triple-grip Ultra Push Chuck enables simple one-touch bur exchange. The auto-lock mechanism automatically controls bur retention when higher torque loads are applied during cutting procedures.

## 4 NEW Ergonomic Arc-Line Body


The distinctive Arc-Line Body constantly offers ideal head and grip angle in any treatment location. While minimising the interruption of your vision, the curved bodyline allows you more freedom of movement and reduces the danger of hitting anterior teeth during complex or delicate clinical procedure.

## Lasting Power & Performance

PANA-MAX is the latest of NSK's popular PANA series.

Every feature of the new PANA-MAX is engineered to deliver unparalleled performance vs. cost ratio, thus helping you work more efficiently and comfortably.

- The highest torque in its class
- Ergonomic Arc-Line Body
- Substantially improved durability
- Lower noise level
- Clean Head System
- NSK Ultra Push Chuck
- Autoclavable at 135°C max
- Thermodisinfectant Ready

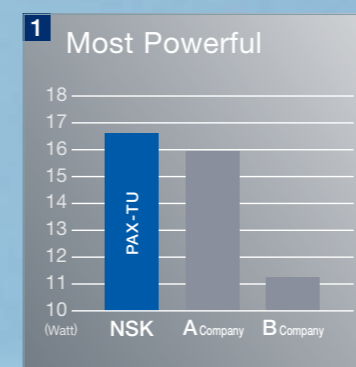

Figures by NSK Laboratory

\*: as of May, 2006

