








SCALING

MODEL	PRODUCTS	ORDER CODE
S1 Universal Type		Z252-411
S2 Sickle Type		Z252-413
S3 Perio Type		Z252-412

PERIODONTICS







S10 Universal Type Slim Line		Z252-420
S20 Root Planing (Straight tip)		Z252-441
S21R Root Planing (Right curved tip)		Z252-442
S21L Root Planing (Left curved tip)		Z252-443




ENDODONTICS V-Tip




S75 V-Tip holder		Z252-035
V-U15 V-file irrigation 23mm #15 ISO15		Y900-177 (Pack of 6)
V-U20 V-file irrigation 23mm #20 ISO20		Y900-178 (Pack of 6)
V-U25 V-file irrigation 21mm #25 ISO25		Y900-179 (Pack of 6)
V-U30 V-file irrigation 22mm #30 ISO30		Y900-180 (Pack of 6)
V-U35 V-file irrigation 23mm #35 ISO35		Y900-181 (Pack of 6)
V-N50 Ni-Ti irrigation 17mm		Y900-182 (Pack of 3)
V-N51 Ni-Ti irrigation 21mm		Y900-183 (Pack of 3)

*Note : Endodontic V-Tips must be used with the NSK S950 series.

RESTORATIVE V-Tip

MODEL	PRODUCTS	ORDER CODE
S75 V-Tip Holder		Z252-035
V-G70 Diamond tip	 φ0.8 (Pack of 3)	Y900-118
V-G71 Diamond tip	 φ1.2 (Pack of 3)	Y900-119
V-G72 Diamond tip	 φ1.6 (Pack of 3)	Y900-120
V-G75 Diamond tip	 φ1.0 φ0.8 (Pack of 3)	Y900-121
V-G76 Diamond tip	 φ1.5 φ0.9 (Pack of 3)	Y900-122
Assorted Tip Kit	includes 5 tips (V-G70, G71, G72, G75, G76)	Y900-158

V-G77 Diamond Tip / Normal Grade	 4.5 φ1.2 (Pack of 3)	Y900-165
V-G78 Diamond Tip / Fine Grade	 4.5 φ1.2 (Pack of 3)	Y900-166
V-G79 Diamond Tip / Super Fine Grade	 4.5 φ1.2 (Pack of 3)	Y900-167
Assorted Tip Kit	includes 3 tips (V-G77, G78, G79)	Y900-168

V-S1 Stainless Tip	 0.9 φ0.8 (Pack of 3)	Y900-161
V-S3 Stainless Tip	 1.2 φ1.1 (Pack of 3)	Y900-162
V-S33 Stainless Tip	 2.6 φ0.8 (Pack of 3)	Y900-163
Assorted Tip Kit	includes 3 tips (V-S1, V-S3, V-S33)	Y900-164

versatile & precise

Designed for a Wide Range of Applications



135°C
Autoclavable at 135°C max

NAKANISHI INC.
www.nsk-inc.com
700 Shimohinata Kanuma-shi
Tochigi 322-8666,
Japan

TEL : +81 (0)289-64-3380
FAX : +81 (0)289-62-5636

NSK Europe GmbH
www.nsk-europe.de
Westerbachstraße 58
D-60489 Frankfurt,
Germany

TEL : +49 (0)69 74 22 99 0
FAX : +49 (0)69 74 22 99 29

NSK United Kingdom Ltd.
www.nsk-uk.com
Office 5, Gateway1000,
Arlington Business Park,
Whittle Way, Stevenage, SG1 2FP, UK

TEL : +44 (0) 1438 310670
FAX : +44 (0) 1438 310671

NSK Oceania Pty Ltd.
www.nsk-inc.com
Unit 22, 198-222 Young St.
Waterloo, Sydney,
NSW 2017 Australia

TEL : +61 (0)2 8306 3000
FAX : +61 (0)2 8306 3030

NSK America - Canadian Branch
www.nsk-inc.com
1054 Centre St. - Suite 283
Thornhill, Ontario, Canada
L4J 8E5

TEL/FAX : +1 905 780 6757
Toll Free : +1 866 780 6757

NSK Middle East
www.nsk-inc.com
Dubai Airport Free Zone,
Room 2W-310,
3F West Wing, Dubai, UAE

TEL : +971 (0)4-299-6644
FAX : +971 (0)4-299-6664



Scaling



Perio



Minimal Intervention



Endodontics



Powerful Partners®


The Power Ring Allows Specific Adjustment for a Variety of Clinical Applications

Scaling / Endodontics / Periodontics / Restorative



- Three level Power Ring
- Solid Titanium body
- Circular Optic Illumination
- Extensive tip range

APPLICATIONS	TIPS	NSK POWER RING SETTING		
		LEVEL 1	LEVEL 2	LEVEL 3
Scaling	S1, S2, S3	←————→		
Periodontics	S10, S20, S21R, S21L	←————→		
Endodontics	S75 + V tips (U-file, Ni-Ti)	←————→		
Restorative	Minimal Intervention	←————→		
	Finishing, Trimming, Polishing	←————→		
	For Caries of Dentin	←————→		

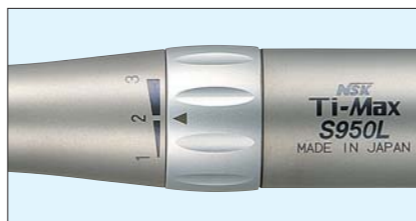


Specific Power Adjustment

NSK Power Ring offers three different power level settings for diversified clinical needs. While the gentle minimum mode is for delicate procedures, the maximum power mode is good for general scaling.



Optic Illumination

Ti-Max S950L model comes with NSK's Circular Optic Illumination. A luminous target spot helps you to make faster and precise treatment.

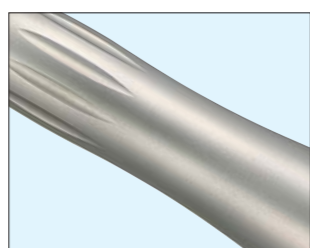

Solid Titanium Body

Crafted from solid titanium with utmost care and precision, S950 is lightweight, ideally balanced and extremely durable. Its titanium body is characterized by a firm and tactile grip, which offers a high degree of operator's control and comfort.



Tip Wrench with Torque Limiter

The tip wrench has a built-in torque limiter device. It eliminates unexpected breakage of tip by preventing over-tightening while exchanging a tip.



The Extensive Tip Range

NSK offers a wide choice of air scaler tips that greatly expand the fields of air scaler applications, including perio, endodontics, scaling and minimal intervention procedures.

Ti-Max
TITANIUM

S950 Air Scalers

Connect to NSK Couplings

NSK PTL-CL-4HV-T 

Ti-Max S950 is a multi-purpose air scaler. In combination with an extensive lineup of NSK air scaler tips, it offers a wide variety of clinical applications. NSK Power Ring allows you to select 3 power levels depending on procedure types such as perio, endodontics, scaling and minimal intervention. With an ergonomically designed body made of solid titanium, S950 is lightweight and extremely durable. Besides NSK couplings, S950 offers convenient direct connectivity to other brands including KaVo®, Sirona®, W&H® and Bien-Air®.

SCALING
PERIODONTICS
ENDODONTICS
RESTORATIVE

Features

Solid Titanium body / Power Ring for 3 power levels / Circular Optic Illumination (S950L) / Frequency of oscillation : 6,200~6,450 Hz



MODEL ORDER CODE

S950L Optic T750

- Weight : 60g
- includes 3 tips (S1, S2, S3), Tip Wrench and Tip Cover

S950 Non-Optic T749

- Weight : 63g
- includes 3 tips (S1, S2, S3), Tip Wrench and Tip Cover

Ti-Max S950KL / S950SL / S950WL / S950BL

Connect Direct to KaVo®, Sirona®, W&H®, Bien-Air® Couplings

SCALING
PERIODONTICS
ENDODONTICS
RESTORATIVE

Features

Solid Titanium body / Power Ring for 3 power levels / Circular Optic Illumination / Frequency of oscillation : 6,200~6,450 Hz



KaVo® MULTIflexLUX® 

S950KL Optic T771

- Weight : 80g
- includes 3 tips (S1, S2, S3), Tip Wrench and Tip Cover



Sirona® Quick Coupling 

S950SL Optic T785

- Weight : 76g
- includes 3 tips (S1, S2, S3), Tip Wrench and Tip Cover



W&H® Roto Quick® 

S950WL Optic T786

- Weight : 84g
- includes 3 tips (S1, S2, S3), Tip Wrench and Tip Cover



Bien-Air® Unifix® 

S950BL Optic T784

- Weight : 85g
- includes 3 tips (S1, S2, S3), Tip Wrench and Tip Cover

*KaVo® and MULTIflexLUX® are registered trademarks of Kaltenbach & Voigt GmbH & Co., Germany. Sirona® is a registered trademark of Sirona Dental Systems GmbH, Germany. W&H® and Roto Quick® are registered trademarks of W&H Dentalwerk Bürmoos GmbH, Austria. Bien-Air® and Unifix® are registered trademarks of Bien-Air Dental S.A., Switzerland.